630-015

DK	Blodtryksmåler til overarm.....	2
SE	Blodtrycksmätare för överarm.....	8
NO	Blodtrykksmåler for overarm.....	13
FI	Verenpainemittari olkavarsimittaukseen.....	18
UK	Blood pressure meter for upper arm.....	23
DE	Oberarm-Blutdruckmessgerät.....	29
PL	Ciśnieniomierz naramienny.....	35

www.adexi.eu

INTRODUKTION

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UK

INTRODUCTION

To get the best out of your new appliance, please read this user guide carefully before using it for the first time. Take particular note of the safety precautions. We also recommend that you keep the instructions for future reference.

IMPORTANT SAFETY PRECAUTIONS

General

- Incorrect use of this appliance may cause personal injury and damage.
- Use for its intended purpose only. The manufacturer is not responsible for any injury or damage resulting from incorrect use or handling (see also Guarantee Terms).
- For domestic use only. Not for outdoor or commercial use.
- Do not immerse the appliance in water and make sure no water enters the appliance.
- Do not expose to direct sunlight, high temperatures, humidity, noise or corrosive substances.
- Never leave unattended when in use, and keep an eye on children.
- Only use accessories manufactured or recommended by the manufacturer.

Use

- This appliance must not be used to make medical diagnoses. It is suitable only for reference. Consult a doctor if your blood pressure readings seem abnormal.
- Always follow the procedure specified in these user instructions when taking readings.
- Do not stretch, twist or bend the cuff or the air hoses.
- Do not use the appliance on individuals who suffer from severe arrhythmia or other cardiac diseases.
- Do not use on babies or others who cannot communicate if the cuff is too tight, etc.
- Do not use in the vicinity of mobile phones or microwave ovens as these emit electro-magnetic radiation that can affect the appliance.
- Do not use if not functioning properly, or if damaged by a heavy impact or similar.


- If the appliance, air hoses or cuff are damaged, have the appliance inspected and if necessary repaired by an authorised repairer. Never try to repair the appliance yourself. Please contact the store where you bought the appliance for repairs under guarantee. Unauthorised repairs or modifications to the appliance will invalidate the guarantee.

Batteries

- Check that the polarities of the batteries are facing in the right direction, as incorrectly fitted batteries can damage the appliance.
- Remove the batteries from the appliance if not used for an extended period of time.
- Always replace all batteries at the same time.
- Never mix different types of battery (e.g. manganese dioxide batteries and alkaline batteries), old and new batteries or different makes of batteries.


KEY TO MAIN COMPONENTS OF THE APPLIANCE

1. End of cuff
2. Cuff
3. Metal bracket
4. Velcro tab
5. Rubber bulb
6. Socket for air hose (not visible on the illustration)
7. Display
8. On/Off button
9. Release valve
10. Hose connector
11. Air hoses
12. Storage pouch

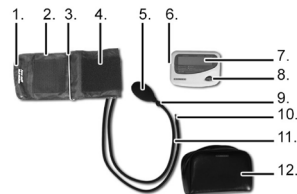Explanation of display and symbols

- a. Systolic pressure
- b. Diastolic pressure
- c. Pulse
- d. Symbol for new measurement
- e. Symbol for pulse measurement
- f. Symbol for emptying of cuff
- g. Battery symbol

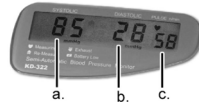d.

e.

f.

ABOUT BLOOD PRESSURE

What is blood pressure?

When the heart pumps blood around the body, the blood exerts pressure on the walls of the blood vessels. This is known as blood pressure.

Blood pressure is different in different parts of the body and depends on the pulse, how elastic the blood vessels are and how thick the blood is.

The pressure that can be measured when the heart pumps out into the body is known as systolic blood pressure. The pressure that can be measured when blood runs back from the body towards the heart is known as diastolic blood pressure.

The two types of blood pressure are usually given together, with the systolic blood pressure always given first. If the blood pressure is shown as "120/70" (pronounced "120 over 70"), it means that the systolic blood pressure is 120 and the diastolic pressure is 70.

The reading unit for blood pressure is mmHg (millimetres of mercury).

Why does blood pressure vary?

Blood pressure varies throughout the day and over a person's life. Many different factors affect your blood pressure, something you should be aware of when taking your blood pressure. For instance, be aware of the following:

- Blood pressure is higher during the day than at night.
- Blood pressure increases slightly after you have just eaten. A lot of salt in food can also increase blood pressure.

- Light exercise (e.g. walking and jogging) causes the systolic blood pressure to rise slightly, while hard exercise will cause it to rise dramatically. However, it will fall again when the body relaxes. A lack of exercise raises blood pressure permanently.
- Blood pressure can increase if you are excited, stressed or tense.
- Sudden climatic changes can cause blood pressure to rise.
- The systolic pressure increases with age.
- Women have a greater risk of increased blood pressure after the menopause.
- Insufficient sleep, smoking, obesity, alcohol and certain illnesses cause blood pressure to rise.
- Bathing and toilet visits may cause blood pressure to rise.
- If the cuff is too loose or too tight or not level with the heart, the measured blood pressure may be lower than the actual blood pressure.

When is my blood pressure too high or too low?

High blood pressure is also called hypertension.

Low blood pressure is also called hypotension.

The table below provides general guidelines for when your blood pressure is too high or too low. You should always consult your doctor for exact readings.

Range	Systolic	Diastolic
Very high blood pressure (hypertension)	> 140 mmHg	> 90 mmHg
High normal	130-139 mmHg	85-89 mmHg
Normal	91-129 mmHg	51-84 mmHg
Erittäin matala Very low blood pressure (hypotension)	< 90 mmHg	< 50 mmHg

Why should I measure my blood pressure?

Blood pressure rises naturally with age. But high blood pressure can also be the first sign of illness, being overweight or of the arteries beginning to stiffen, which increases the risk of blood clots, among other things. It is therefore a good idea to keep an eye on your blood pressure and consult a doctor if it seems unusually high or low. You should write down your measurements so you can show them to your doctor.

Important! Do not use the appliance for medical diagnosis – only for reference purposes. Consult a doctor if your blood pressure readings seem abnormal. Always follow the procedure specified in these user instructions when taking readings.

PRIOR TO FIRST USE

Inserting the batteries

- Open the battery cover at the back of the appliance by pushing the cover in the direction of the arrow and removing it.
- Insert four batteries, size “AA”. Make sure the batteries are fitted with the correct polarity (see the markings at the bottom of the battery compartment).
 - Do not use rechargeable 1.2 V batteries for the appliance.
- Replace the cover.
- Replace the batteries when the battery symbol (g) appears on the display (7).

Connecting the air hose

- Connect hose connector (10) to air hose hole (6).
- Check that the hose connector is inserted as far as possible into the appliance and is firmly in place.

USING THE APPARATUS

Preparing for a blood pressure reading

- Sit or lie down for 5-10 minutes before the reading.
- Avoid eating, smoking, exercising, bathing or drinking alcohol, coffee or tea before taking the reading.
- Do not go to the toilet immediately before the reading.
- Do not take a blood pressure reading if you are excited, stressed or tense.

- Relax during the reading and do not talk.

Fitting the cuff on the arm

- Thread the end of the cuff (1) through the metal bracket (3) so the cuff (2) forms a tube. Make sure the Velcro tabs face outwards. Turn the cuff so the air hoses (11) point towards the floor.
- Slide the cuff onto the arm and position it on the upper arm, approx. 1-2 cm above the elbow.
 - The cuff should be approximately level with your heart to obtain an accurate reading.
- Tighten the cuff by pulling at the end until the cuff fits tightly around the arm. At the same time, turn the cuff so the air hoses run along the inside of the arm (see illustration).
 - The hoses should be parallel to your arm and point towards your hand.
 - Do not overtighten the cuff – you should be able to fit one finger between the cuff and your arm.
- Fold the end of the cuff over the metal bracket, wrap the end down around your arm and press it onto the Velcro tab (4).
- NB:
 - The cuff should be taut but should not pinch.
 - The air hoses should run along the inside of your arm.
 - The lower edge of the cuff should be 1-2 cm from your elbow joint.
 - It is best to place the cuff directly against the skin.
 - If you want to take a number of readings, place the cuff on the same arm each time. It is also important to also always measure at the same time of day, under the same conditions and in the same position.

Measuring the blood pressure

- You can take the reading in one of two positions:
 - Lie down on your back with your arm resting along your side and the palm facing upwards.
 - Sit on a chair and place your arm on the table with your palm facing upwards.

Stay completely still while taking the reading.

- Push the On/Off button (8). The appliance beeps and all symbols and numbers are shown on the display for approx. 2 seconds.
- When the display shows “0 mmHg” under both systolic pressure (a) and diastolic pressure (b) and the symbol for new reading (d) flashes, the appliance is ready for use.
 - If air was left in the cuff from last use, the symbol for emptying the cuff (f) flashes. Press the release valve (9) on the hose in front of the rubber bulb (5) to empty all air out of the cuff. When the symbol for new measurement flashes, the appliance is ready for use.
- Inflate the cuff to 30 mmHg above the expected (diastolic) blood pressure by squeezing the rubber bulb a number of times in quick succession. For guidelines about blood pressure values, see the table in “When is my blood pressure too high or too low?” It is important to inflate the cuff in one continuous process without pauses. Let go of the rubber bulb once the desired value has been reached.
 - If an error occurs during inflation of the cuff, the display shows “EE RETEST”. Empty all air out of the cuff by pressing the release valve. Inflate the cuff again when the symbol for new measurement flashes.
- The blood pressure measurement begins shortly after the rubber bulb is released.
 - The display shows “TESTING” and the appliance displays the pressure continuously during the measurement.
 - The symbol for pulse measurement (e) flashes as an indication that the pulse is also being measured.
- When the blood pressure measurement is complete, the appliance beeps three times.
 - The display shows the systolic pressure, the diastolic pressure and the pulse (c).
 - The symbol for emptying of cuff flashes.
- Press the release valve to empty all air out of the cuff. When all air has been removed from the cuff, the symbol for new measurement flashes.
- If you need to take several readings, you should wait 5-10 minutes between each reading.
 - Caution! A single reading will not give a sure indication of your condition. You can only get an idea of this after a number of readings.
- Switch off the appliance by pressing the On/Off button when you have finished the measurement.
 - The appliance automatically switches off after two minutes of inaction.

TROUBLESHOOTING

Error	Cause	Solution
The display shows "RETEST"	<ol style="list-style-type: none"> 1. The pressure in the cuff is not 30 mmHg higher than the expected blood pressure. 2. The user moved the arm or the appliance during measurement. 3. The pressure drops fast due to leakage in the circuit. 	<ol style="list-style-type: none"> 1. Take a new reading and make sure the cuff is inflated to 30 mmHg above the expected blood pressure. 2. The user must remain completely still and not move the arm or the appliance during measurement. 3. Check for leakage in the circuit.
The pressure only rises slowly during inflation.	The hose is incorrectly fitted.	Fit the hose correctly.
The display shows the battery symbol or the display is blank.	The batteries are almost dead.	Replace the batteries. Check that the batteries are fitted the right way round.
Unusual readings	<ol style="list-style-type: none"> 1. The cuff has not been tightened sufficiently or has been positioned incorrectly. 2. The user moves or speaks during measurement. 3. The user is too nervous or tense during measurement. 4. The user is not in the correct lying or sitting position. 5. The user pauses during inflation or uses the appliance incorrectly. 6. The user may have a heart problem. 	<ol style="list-style-type: none"> 1. Tighten the cuff and position it as described in the section "Fitting the cuff on the arm". 2. The user must remain still and avoid speaking during measurement. 3. The user must try to relax and take 2-3 deep breaths. 4. The user must lie down or sit as described in "Measuring the blood pressure". 5. Inflate the cuff as described in "Measuring the blood pressure". 6. Check result by performing one or more new measurements. Talk to your doctor if you consistently have abnormal readings.

STORAGE

Store the appliance in the accompanying storage pouch (12) when not in use.

- Store in a dry, cool place.
- Do not place heavy objects on top of the appliance.
- Remove the batteries if the appliance is not to be used for a prolonged period.
- Do not stretch, twist or bend the air hoses or the cuff.

CLEANING

When cleaning the appliance, you should pay attention to the following points:

- Do not submerge the appliance in water, and do not let water enter the appliance.
- Clean the appliance and the cuff by wiping them with a damp cloth. A little detergent can be added to the water if the appliance and the cuff are heavily soiled.
- Do not use a scouring pad, steel wool or any form of strong solvents or abrasive cleaning agents to clean the appliance, as they may damage the outside surfaces of the appliance.

- **IMPORTANT!** Allow to dry fully after cleaning before using again. Do not use appliance if damp.

INFORMATION ON DISPOSAL AND RECYCLING OF THIS PRODUCT

Please note that this Adexi product is marked with

this symbol: 

This means that this product must not be disposed of along with ordinary household waste, as electrical and electronic waste must be disposed of separately.

In accordance with the WEEE directive, every member state must ensure correct collection, recovery, handling and recycling of electrical and electronic waste. Private households in the EU can take used equipment to special recycling stations free of charge. In certain member states, used apparatus can be returned to the dealer where they were bought on the condition you buy new products. Contact your retailer, distributor or the municipal authorities for further information on what to do with electrical and electronic waste.

GUARANTEE TERMS

The guarantee does not apply:

- if the above instructions are not followed
- if the appliance has been interfered with
- if the appliance has been mishandled, subjected to rough treatment, or has suffered any other form of damage
- if faults have arisen because of faults in your electricity supply.

Due to the constant development of our products in terms of function and design, we reserve the right to make changes to the product without prior warning.

FREQUENTLY ASKED QUESTIONS

If you have any questions regarding use of the appliance and cannot find the answer in this user guide, please try our website at www.adexi.eu. Go to the "Consumer Services" menu, click on "Question & Answer" to view the most frequently asked questions.

You can also see contact details if you need to contact us for technical issues, repairs, accessories and spare parts.

IMPORTER

Adexi Group
www.adexi.eu

We cannot be held responsible for any printing errors.

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